

EZ-ACCESS[®] Product Catalog

Anytime, anywhere you need access



6

Residential Access
Product Lines

Enriching lives
by providing
access to life
beyond barriers

Family-owned and
-operated since

1984



www.ezaccess.com
customerservice@ezaccess.com
1-800-451-1903

Solution-based products to make safe, in



1

PATHWAY® 3G Modular Access System:
Create safe and seamless access for both standard and unique entryways with this semi-permanent, customizable ramp system.

GATEWAY® 3G Solid Surface Portable Ramp:
Perfect for temporary or semi-permanent home access, the GATEWAY 3G can be easily moved or relocated as needs change. Available with or without handrails.


2

TRANSITIONS® Threshold Ramps:
Create a smooth transition and eliminate tripping hazards with this convenient solution.

3

TILT® Toilet Incline Lift:
Experience both comfort and safety with this innovative toilet aid that offers safe, independent use and assistance in getting on and off the toilet.

Independent living a possibility

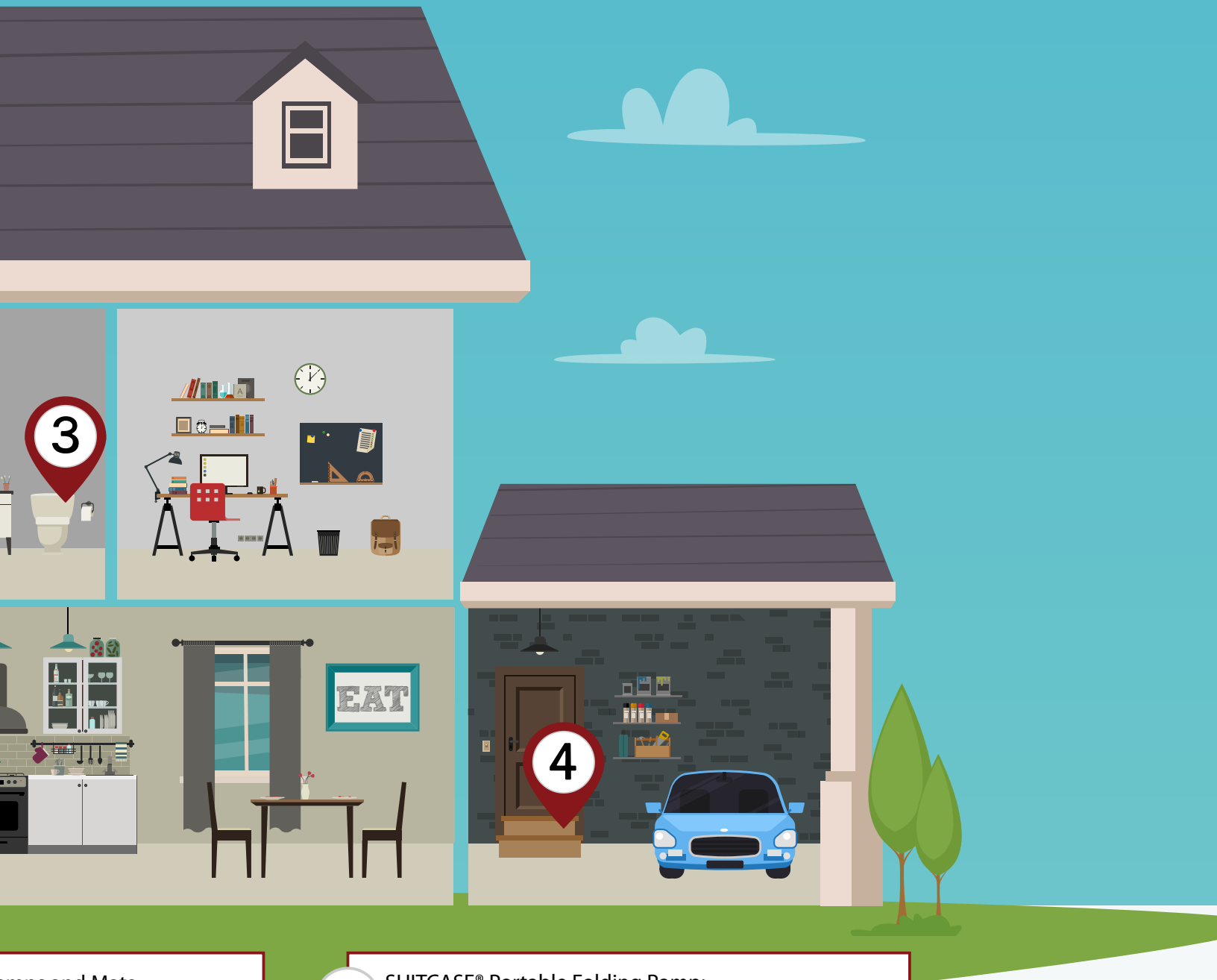

Ramps and Mats:
Prevent falls and eliminate potential trip
hazard solutions.

Improve safety with this toileting
solution for independent bathroom use by
removing the need to get off of the commode.

4

SUITCASE® Portable Folding Ramp:
Overcome curbs and low landings with this versatile
portable ramp that provides the safety and strength
required while remaining lightweight and portable.
Great for at home or on the go.

PASSPORT® Vertical Platform Lift:
The ideal solution when surrounding space is
limited or higher rises are present. Great for use at your
front or back door, as well as garage entry.

SUITCASE® Portable Ramps

SUITCASE® ramps are a durable and dependable solution for a variety of mobility and accessibility needs. From the SUITCASE Singlefold to the Singlefold AS and Trifold AS, you are sure to find a convenient, lightweight, and easy-to-use ramp perfect for wheelchair and scooter use. These portable ramps are built to withstand repeated use and offer a variety of slip-resistant surfaces. All SUITCASE ramps feature non-breakable handles are ideal for use at home or on the go.

SUITCASE® Singlefold Ramp



Weight capacity		800 lbs.
Width		30"
Material		Aluminum
Tread		Extruded
Warranty		Lifetime

Made in the USA

- » Single-fold design sets up quickly
- » Full-length hinge for additional strength and stability
- » Made of lightweight, durable, and high-strength aluminum
- » Extruded, slip-resistant surface for superior traction

SUITCASE® Singlefold AS Ramp



Weight capacity		800 lbs.
Usable width		29.5"
Material		Aluminum
Tread		Applied
Warranty		Lifetime

Made in the USA

- » Single-fold design sets up quickly and folds for storage or transport
- » Can be used as one unit or easily separated into two lightweight sections
- » Made of lightweight, durable, and high-strength aluminum
- » Applied, slip-resistant surface

SUITCASE® TRIFOLD® AS Ramp



5'

6'

7'

8'

10'

Weight capacity		800 lbs.
Usable width		29.5"
Material		Aluminum
Tread		Applied
Warranty		Lifetime

Made in the USA

- » Three-fold hinge allows for ultimate portability when longer ramp lengths are needed
- » Designed to be used and carried as one unit, or separated into two individual sections for easier transport and reduced handling weight
- » Made of lightweight, durable, and high-strength aluminum
- » Applied, slip-resistant surface

SUITCASE® Top Lip Extension



Weight capacity		500 lbs.
Material		Aluminum
Tread		Extruded
Warranty		Lifetime

Increase the reach of SUITCASE ramps with the optional Top Lip Extension. The TLE extends the standard 3" top transition plate to 9", allowing the ramp to clear the rear bumper when loading a mobility device into the back of an SUV or van.

- » Extruded, slip-resistant surface
- » Made of durable, high-strength aluminum
- » Sold as a pair

TRANSITIONS® Threshold Ramps

Create a smooth transition for mobility users and eliminate potential trip hazards with TRANSITIONS® threshold ramps, mats, and plates. This convenient solution is available in a variety of surface and material options - as well as height ranges - and accommodates most threshold applications.

Modular Entry Ramp

- » Modular design with interchangeable extensions
- » Available in seven sizes to accommodate thresholds from 3/4" to 6-1/4" high
- » Made of high-strength aluminum which will not rust or corrode
- » Extruded, slip-resistant surface for superior traction
- » 34" wide
- » 700 lb. weight capacity



Angled Entry Ramp

- » Independently adjustable legs with swivel feet to fit varying heights and angles
- » Three lengths available with vertical adjustments up to 5-7/8" high
- » Made of high-strength aluminum which will not rust or corrode
- » Extruded, slip-resistant surface for superior traction
- » 36" wide
- » 700 lb. weight capacity




Angled Entry Plate

- » Single-piece, formed plate with an applied, slip-resistant surface
- » Available in 8", 10", and 12" lengths
- » Accommodates thresholds up to 2" high
- » Made of high-strength aluminum that will not rust or corrode
- » 32" wide
- » 700 lb. weight capacity



Angled Entry Mat

- » Made of recycled rubber
- » Beveled edges for cross-traffic access
- » Available in 1-1/2" and 2-1/2" heights and can easily be trimmed or notched if needed
- » Offered with or without the recessed "Welcome"
- » Black or gray color options
- » Usable width of 36"
- » 850 lb. weight capacity

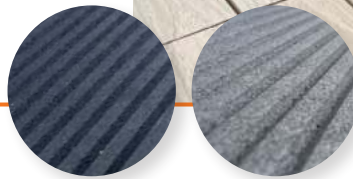

Modular Entry Mat


- » Made of 100% recycled rubber
- » 2-1/2" height or combine with optional risers to obtain a maximum height of 4-3/4"
- » Can be trimmed or notched for a custom fit
- » Ideal for both indoor or and outdoor use
- » 48" wide
- » 850 lb. weight capacity



All TRANSITIONS products are backed by our lifetime warranty.

PATHWAY® 3G Modular Access System

The American-made PATHWAY® 3G Modular Access System consists of ramps, stairs, and platforms, all of which may be combined and configured in a variety of ways to provide a customized home access solution. The ramp's all-aluminum, modular design offers benefits to both users and installers and has the ability to create a safe pathway for individuals of all mobility levels.

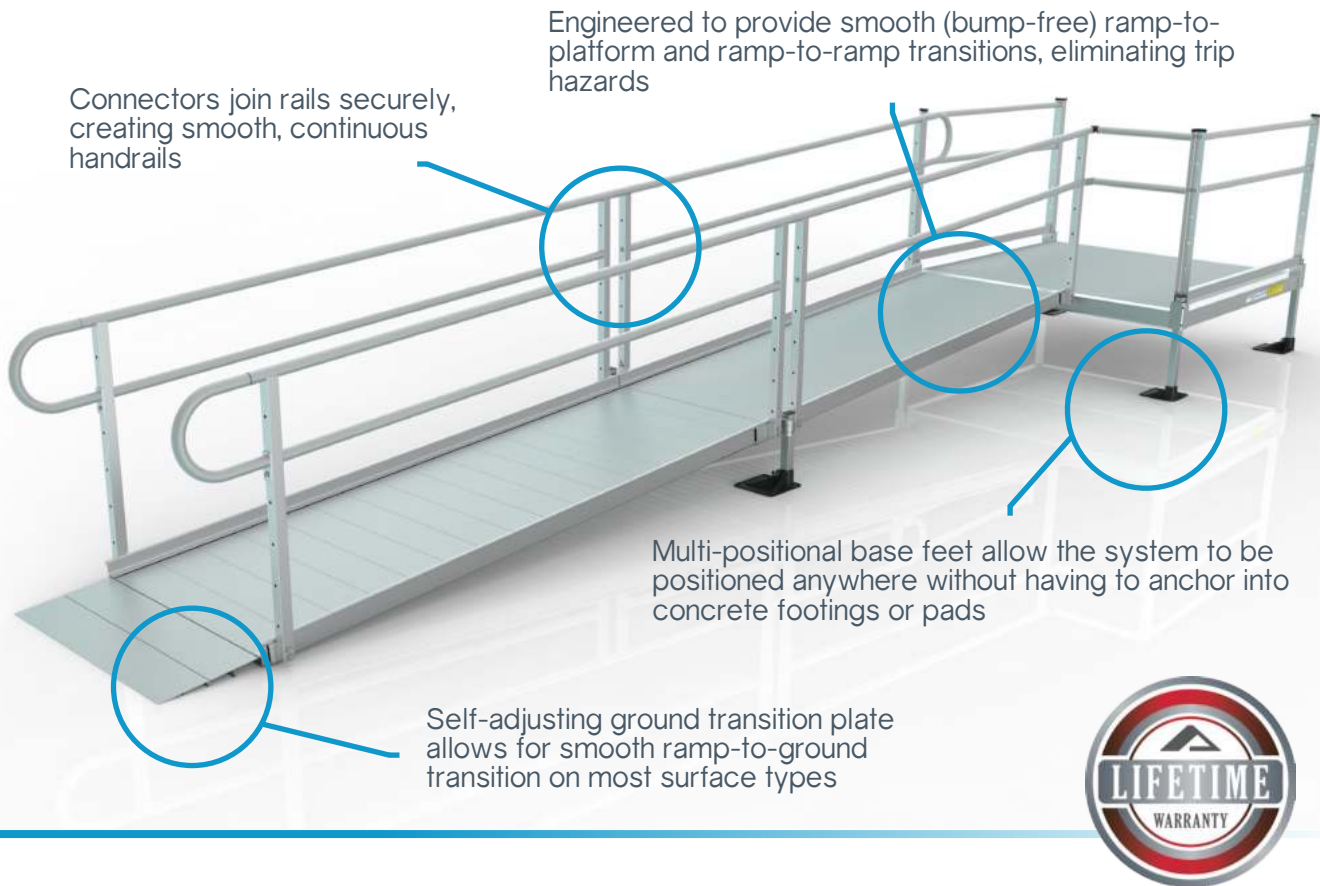

- » Made of lightweight, durable, and high-strength aluminum that will not rust, making this system virtually maintenance-free
- » Can be easily reconfigured or relocated as needs change
- » Quick installation with minimal tools
- » Interchangeable surface and handrail options for immense versatility
- » No poured footers, no construction - just simple assembly

Ramps



- » Weight capacity: 1,000 lbs.
- » Standard widths: 36" & 48" (48" width available in solid surface only)
- » Ramp handrail height: 34-1/2"


Platforms




- » Weight capacity: 1,000 lbs.
- » Platform handrail height: 34-1/2"

Surfaces

Surface options allow for infinite ramp configurations and design flexibility while providing safe access to every home. Both surfaces are made of durable aluminum and offer advanced slip-resistance.

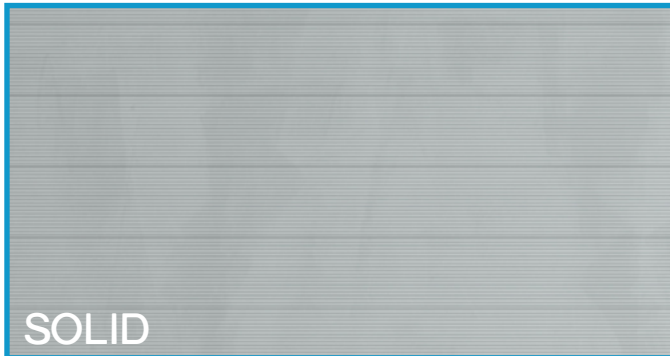

The solid surface, our most popular option, has an extruded tread and features raised ribs to provide slip-resistance and all-weather traction.

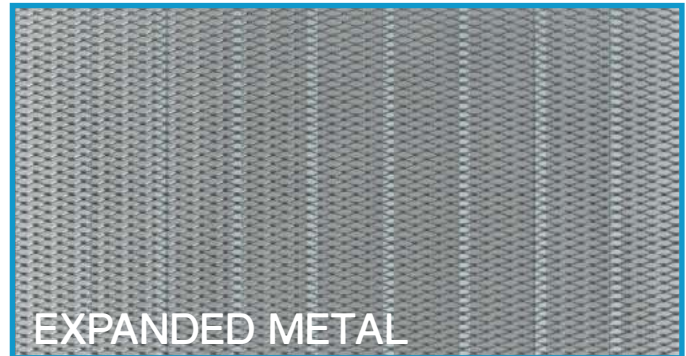

The diamond-shaped openings in the expanded metal surface allows multi-directional traction and reduced surface accumulation by allowing small debris, rain, and melting snow to flow through. Only available in select ramp and platform sizes.

Handrails

A variety of handrail styles are available to suit your need. Whether you choose a two-line rail, vertical pickets, or horizontal pickets, these versatile, component-based rails are easy to assemble and customize. All styles are made of aluminum and feature a durable design with continuous top rails.



» **Two-line Rail:** The most popular handrail design, the two-line rail provides an additional gripping surface by adding a second rail below the continuous top rail.



» **Vertical Pickets:** Vertical pickets offer an added level of safety, especially for children and small pets. Our pre-assembled aluminum inserts meet the 4-inch sphere rule for spacing.



» **Horizontal Pickets:** Horizontal pickets offer an aesthetic option while maintaining safety. These horizontal inserts run parallel to the topmost rail and are only available with solid surface ramps.

GATEWAY® 3G Solid Surface Portable Ramp

The GATEWAY® 3G Solid Surface Portable Ramp is the ideal solution when mobility needs extend beyond the scope of a portable folding ramp. Designed for the rigors of repeated, everyday use in all weather conditions and is ideal for scooters, wheelchairs, and walkers. Available with or without handrails.



- 3'
- 4'
- 5'
- 6'
- 7'
- 8'
- 9'
- 10'



- » More durable than wood, the GATEWAY 3G is made of aluminum and designed for years of maintenance-free use
- » Solid, seam-free surface that offers slip-resistance
- » Perfect for temporary or semi-permanent home access and can be easily moved or relocated as needs change
- » Formed top transition for secure placement
- » Adjustable bottom transition that allows smooth ramp-to-ground conversion
- » 36" wide with raised side rails to help prevent mobility devices from driving off the edge
- » 1,000 lb. weight capacity

RAMP ONLY

RAMP WITH TWO-LINE HANDRAILS

RAMP WITH VERTICAL PICKET HANDRAILS



The GATEWAY 3G is backed by our lifetime warranty and is made in the USA.

PASSPORT® Vertical Platform Lift

The PASSPORT® Vertical Platform Lift takes safe access to the next level, providing a quiet, smooth ride to your next destination. When surrounding space is limited or higher rises are present, the PASSPORT is the perfect solution. The PASSPORT is user-friendly, low-maintenance, and is the preferred lift of technicians and installers. ETL certified – independently tested and approved for product safety.



Highlighted Features

- » 52" lifting height available with straight or turn platforms
- » 750 lb. weight capacity
- » Lightest residential platform lift on the market
- » Plug-N-Play wiring technology for quick and easy installation
- » DC-powered unit
- » On-board smart charger plugs into standard GFI outlet and offers approximately 30 complete, fully loaded cycles on battery power alone
- » High-visibility LED diagnostic indicator lights for troubleshooting and reduced service calls
- » Multi-positional control box features keyed operation
- » Continuous, large, low-pressure mushroom head control switches (including emergency stop) are easy to use, regardless of motor skills
- » Undercarriage safety pan with obstruction detector (takes less than 15 pounds of pressure to stop lift)
- » Self-actuating, 24-inch entry ramp folds upward when the lift rises to prevent accidental travel off the platform

The PASSPORT is proudly made in the USA and backed by our 2-year warranty.

TILT® Toilet Incline Lift

The TILT® Toilet Incline Lift is a patented toileting aid that lifts to your natural body movement, providing the independence needed to maintain dignity in the bathroom. By combining the functionality of a heavy duty commode and a lift chair, the lightweight and durable TILT offers an easy-to-use bathroom safety solution. It also features handheld control to make the lifting process convenient for the user or caregiver. Designed for quick, simple installation. Available in corded or battery power.

Corded Power Unit



Battery Powered Unit



- » Lightweight and features a 325 lb. weight capacity
- » Can be installed in as little as 15 minutes
- » Standard or elongated seat models available



Equipped with a Companion Control to allow the user or caregiver to easily operate the seat with the push of a button.



Equipped with non-slip hand grips, which allows users to gain a secure grip while the seat raises and lowers.



Battery powered unit features a rechargeable battery pack that mounts to the wall within 3' of the unit and requires 5" x 16" of wall space.



The TILT bolts directly to the bowl for increased stability.

The TILT is backed by our 2-year warranty.





www.ezaccess.com
customerservice@ezaccess.com
1-800-451-1903